



Genotyping protocol

Vat1l

IR00002809 / E127

(ICS internal reference)

This report has been prepared by: **Christelle Roth**
genotyping@igbmc.fr

This report has been validated by: **Sylvie Jacquot, PhD, Head of Genotyping Service**
33 (0)3 88 65 57 44
genotyping@igbmc.fr

The first version of this report was generated the: 19 May 2014

For any question, please contact:

Institut Clinique de la Souris - ICS - Mouse Clinical Institute
1 rue Laurent Fries, BP 10142
67404 Illkirch Cedex, France
Email: genotyping@igbmc.fr
Web site: <http://www-mci.u-strasbg.fr/>

TABLE OF CONTENTS

Table of contents2

1. Genotyping protocol and data2

 1.1. Genotyping strategy2

 1.2. PCR protocol4

2. Cre and Flp genotyping method5

1. Genotyping protocol and data


This section describes the condition used at the Mouse Clinical Institute (ICS) to genotype your **Vat1l** Constitutive Knockout / Conditional Knockout (KO-cKO) project.

1.1. Genotyping strategy

The map below describes the position of the primers used for genotyping for each possible allele.



KO-cKO pl _ Genotyping strategy

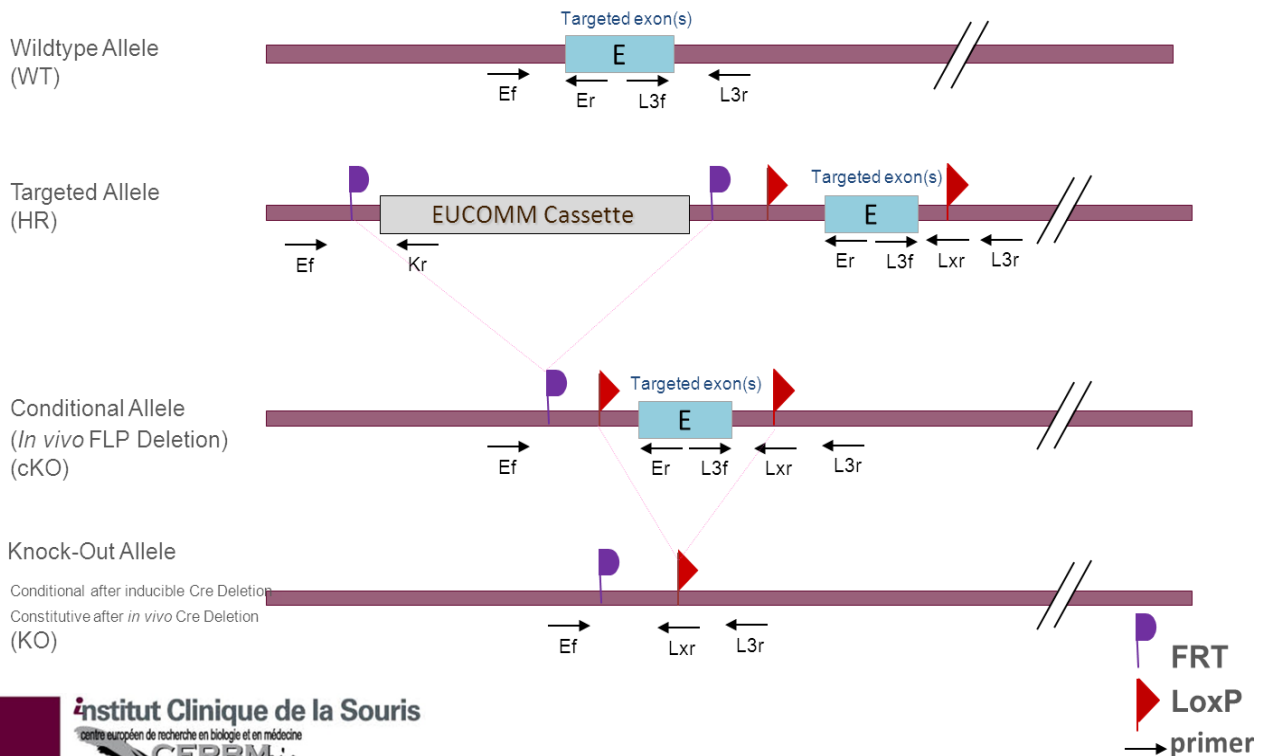

Sequence of primers used for genotyping:

Position	Primers	Sequence
Ef	5096	GATGGCTCGCGTTCATTAT
Er	5100	GACATCCCAGCACTTAATCAAAGGCTT
Kr	3277	CTCCTACATAGTTGGCAGTGTTTGGG
L3f	5098	AGGATTCGAATGTTCTGGGA
L3f ²	5099	TTCTAGAGCCTGTGAGGCAGTAAAAGT
L3r	5097	CCTGCTCTATCCATTCTAAAATA
Lxr	3255	ACTGATGGCGAGCTCAGACCATAAC

²: for a selected position, a second primer was designed

PCR fragments expected size (bp):

Region analyzed	Primers used	Position on the primer (see the map above)	Targeted allele (HR)	conditional allele (KO-cKO)	KO allele	WildType allele
5' part of the selection marker	5096-3277	Ef / Kr	314	---	---	---
Presence of the distal loxP (with DMSO)	5098-5097 (with 5% DMSO)	L3f / L3r	384	384	---	384
Distal loxP specific PCR	5099-3255	L3f ² / Lxr	237	237	---	---
Excision of the selection marker	5096-5100	Ef / Er	7310*	406	---	227

*: this PCR product will not be observed using our PCR genotyping conditions (see description below)

** : this PCR is only verified if mice are generated

---: no Amplicon should be obtained

1.2. PCR protocol

This section describes the composition of the mix and cycling conditions used for genotyping.

Reagents:	Volume:
- FastStart PCR Master (Roche)	7.5µl
- DNA (50ng/µl)	1.5µl
- 5' primer (100 µM)	0.06µl
- 3' primer (100 µM)	0.06µl
- Sterile H ₂ O	up to 15 µl

Cycling conditions:

Temp	Time	#Cycles
95°C	4min	1
94°C	30s	34
62°C	30s	
72°C	1min	
72°C	7min	1
20°C	5min	1

NB: These PCR conditions have been optimized for high-throughput genotyping. Adaptation to small-scale may be required.

2. Cre and Flp genotyping method

You will find the genotyping protocol in the publication on:

[Highly-efficient, fluorescent, locus directed cre and FlpO deleter mice on a pure C57BL/6N genetic background.](#)

Birling MC, Dierich A, Jacquot S, Hérault Y, Pavlovic G.
Genesis. 2012 Jun;50(6):482-9. doi: 10.1002/dvg.20826. Epub 2012 Mar 20.