

Gene: Zup1


Colony prefix: DAKJ

Allele: Zup1^{em1(IMPC)Wtsi}

Allele type: Crispr/Cas9 mediated deletion

Allele information:

Further information about the allele can be found on the 'International Mouse Phenotyping Consortium' (IMPC) web site at http://www.mousephenotype.org/data/alleles/MGI:1919830/em1%2528IMPC%2529Wtsi

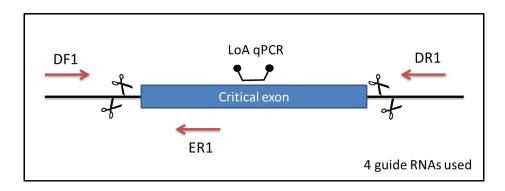

Mouse QC information

Loss of WT Allele (LOA) qPCR	Pass	Mutation Sequence confirmed	Pass
Mutant Specific SR-PCR	Pass	Off-target analysis complete	na

Mutant Allele sequence:

Deletion size (bp): 500

This technical data sheet and information ("Datasheet") is supplied by Genome Research Limited ("GRL").

Although reasonable care is taken in the preparation of this Datasheet, GRL gives no warranties express or implied for any use of the Datasheet or for the accuracy of the Datasheet. GRL assumes no responsibility or liability for any decisions based upon the Datasheet. Without limiting the foregoing the Datasheet was prepared for mice supplied directly from GRL and where copies of this Datasheet are available from third party repositories or distribution centres ("Third Parties") GRL shall not be liable for any inconsistency between the mouse strain supplied by the Third Party and the Datasheet howsoever arising.

Report Generated on: 11-04-2018

Wellcome Sanger Institute
Wellcome Trust Genome Campus

MGPgenotyping@sanger.ac.uk

Hinxton, Cambridge CB10 1SA, UK

www.sanger.ac.uk



Guide RNAs used in initial experiment

Wellcome Sanger Institute Wellcome Trust Genome Campus Hinxton, Cambridge CB10 1SA, UK

MGPgenotyping@sanger.ac.uk www.sanger.ac.uk

Sequence	Chr	Chr Start	Chr End
CAAGCAAGAATGGGATTGCCTGG	10	33939808	33939830
GTTATAAGCTGCCTGAAGTAGGG	10	33939848	33939870
CCAGTCAGTGCTTATGCGCTAAC	10	33940304	33940326
CCTTAACTAGCGGATGTGAACAT	10	33940337	33940359

Genotyping by end-point PCR

This technical data sheet and information ("Datasheet") is supplied by Genome Research Limited ("GRL").



Wellcome Sanger Institute Wellcome Trust Genome Campus Hinxton, Cambridge CB10 1SA, UK

MGPgenotyping@sanger.ac.uk www.sanger.ac.uk

These mice may be genotyped through a combination of separate PCR reactions that detect the gene-specific wild type allele and a mutant allele-specific short

range PCR. Interpretation of the consolidated results produces the genotype of the mice. In addition to the expected product, the mutant assay may also amplify the endogenous wild type sequence which will appear as a larger band on an agarose gel. The presence of this extra band will depend on the size of the original deletion.

PCRs primer pairs and expected size bands

Assay Type	Assay	Forward Primer	Reverse Primer	Expected Size Band (bp)
Standard PCR	Wild type	Zufsp_DF1	Zufsp_ER1	383
Standard PCR	Mutant	Zufsp_DF1	Zufsp_DR1	500

Primer sequences

Primer Name	Primer Sequence (5' > 3')
Zufsp_DF1	CCTGTGGGTTTCCTGTGAAG
Zufsp_ER1	GGCAAAGGGCCTAGATTCAT
Zufsp_DR1	AAACAGTTGGAGGTTGGGAAT

Reaction setup

Reagent	μl
DNA (~50-100 ng)	1
10x Buffer	1.5
MgCl2 (50 mM)	0.45
Platinum Taq (Invitrogen)	0.15
dNTPs (100 mM)	0.15
Primer 1 (10 µM)	0.3
Primer 2 (10 µM)	0.3
ddH20	11.15
Total	15

Amplification conditions

Step	Conditions	Time
1	94°C	5 min
2	94°C	30 sec
3	58°C	30 sec
4	72°C	1:30 sec
5	Go to '2' + 34 cycles	-
6	72°C	5 min
7	12°C	forever

This technical data sheet and information ("Datasheet") is supplied by Genome Research Limited ("GRL").



Genotyping by loss of WT allele qPCR Assay (gene-specific assay)

Wellcome Sanger Institute Wellcome Trust Genome Campus Hinxton, Cambridge CB10 1SA, UK

MGPgenotyping@sanger.ac.uk www.sanger.ac.uk

The wild type loss of allele (LoA) qPCR assay uses a hydrolysis probe assay (for example Applied Biosystems TaqMan® technology) to determine the copy number of the wild type allele in a sample. Homozygotes will show no amplification, heterozygotes one copy and wild type mice will show two copies when compared to a wild type control.

The number of copies of the wildtype allele can be detected using a FAM-labelled custom qPCR TaqMan® assay. These are multiplexed with a VIC® labelled endogenous control assay (for example TaqMan® Copy Number Reference Assay, Mouse, Tfrc; Applied Biosystems part #4458366). Reference DNA controls of known genotypes should also be included to facilitate correct analysis.

Primers for LoA qPCR assay

Gene	Source	Forward Primer Seq.	Reverse Primer Seq.	Probe Primer Seq.
Zufsp	Life Technologies	AAAGCTAAGGTTCACAGTCA AAACAA	CCCCAATAAAATATAACGTACACA GATAA	TATAACATGGACAGACTAGTAAA

Reactions are performed in a 10µl volume using an Applied Biosystems 7900HT Fast Real-Time PCR System or Applied Biosystems Viia7 with DNA prepared using the Sample-to-SNPTM kit (Applied Biosystems) from mouse ear biopsies. GTXpressTM buffer is also used (Applied Biosystems).

Reagent	μl
2x GTXpressTM buffer	5
20x target assay	0.5
ddH2O	3
Tfrc endogenous 20x assay	0.5
DNA	1

Amplification conditions

Step	Conditions	Time
1	95°C	20 sec
2	95°C	10 sec
3	60°C	30 sec
4	Go to '2' + 34 cycles	-

Links to information and frequently asked questions

This technical data sheet and information ("Datasheet") is supplied by Genome Research Limited ("GRL").



Hinxton, Cambridge CB10 1SA, UK

Wellcome Sanger Institute
Wellcome Trust Genome Campus

MGPgenotyping@sanger.ac.uk www.sanger.ac.uk

MGP mouse phenotype data: http://www.mousephenotype.org

Useful publications

White, J.K., Gerdin, A.-K., Karp, N.A., Ryder, E., Buljan, M., Bussell, J.N., Salisbury, J., Clare, S., Ingham, N.J., Podrini, C., et al. (2013). Genome-wide Generation and Systematic Phenotyping of Knockout Mice Reveals New Roles for Many Genes. Cell 154, 452–464.

Mali P, Yang L, Esvelt KM, et al (2013) RNA-guided human genome engineering via Cas9. Science 339:823–6. doi: 10.1126/science.1232033

Jinek M, Chylinski K, Fonfara I, et al (2012) A programmable dual-RNA-guided DNA endonuclease in adaptive bacterial immunity. Science 337:816–21. doi: 10.1126/science.1225829

Cong L, Ran FA, Cox D, et al (2013) Multiplex genome engineering using CRISPR/Cas systems. Science 339:819–23. doi: 10.1126/science.1231143

Singh P, Schimenti JC, Bolcun-Filas E (2014) A Mouse Geneticist's Practical Guide to CRISPR Applications. Genetics genetics.114.169771–. doi: 10.1534/genetics.114.169771

Brandl C, Ortiz O, Röttig B, et al (2015) Creation of targeted genomic deletions using TALEN or CRISPR/Cas nuclease pairs in one-cell mouse embryos. FEBS Open Bio 5:26–35. doi: 10.1016/j.fob.2014.11.009

Zhou J, Wang J, Shen B, et al (2014) Dual sgRNAs facilitate CRISPR/Cas9 mediated mouse genome targeting. FEBS J. doi: 10.1111/febs.12735

Kraft K, Geuer S, Will AJ, et al (2015) Deletions, Inversions, Duplications: Engineering of Structural Variants using CRISPR/Cas in Mice. Cell Rep. doi: 10.1016/j.celrep.2015.01.016

Shen B, Zhang J, Wu H, et al (2013) Generation of gene-modified mice via Cas9/RNA-mediated gene targeting. Cell Res 23:720–3. doi: 10.1038/cr.2013.46

Wang H, Yang H, Shivalila CS, et al (2013) One-step generation of mice carrying mutations in multiple genes by CRISPR/Cas-mediated genome engineering. Cell 153:910–8. doi: 10.1016/j.cell.2013.04.025

Yang H, Wang H, Shivalila CS, et al (2013) One-Step Generation of Mice Carrying Reporter and Conditional Alleles by CRISPR/Cas-Mediated Genome Engineering. Cell 154:1370–1379. doi: 10.1016/j.cell.2013.08.022

This technical data sheet and information ("Datasheet") is supplied by Genome Research Limited ("GRL").